



ONE PLANET AESTHETICS

☎ 0800 772 0402

🌐 www.oneplanetaesthetics.co.uk

✉ info@oneplanetaesthetics.co.uk

📷 📘 @oneplanetaesthetics



St Clare House, Ninth Floor, Greyfriars,
Princes Street, Ipswich, Suffolk, IP1 1LP



ELEVATING AESTHETICS WITH PURPOSE, PRECISION & INTEGRITY

At One Planet Partnership, we are more than a technology distributor — we are passionate curators of transformative skin and body solutions. Founded on over three decades of industry experience, our mission is clear: to connect clinics, aesthetic professionals, and medical practitioners with cutting-edge devices that are effective, ethical, and future-focused.

OUR CORE VALUES ARE ROOTED IN:

Clinical Integrity: We only represent technologies that deliver measurable, evidence-based results.


Partnership with Purpose: We build long-term relationships with professionals who share our commitment to safety, results, and client satisfaction.

Innovation with Intention: Every device in our portfolio is handpicked for its quality, innovation, and real-world performance.

We believe in creating a more intelligent and responsible aesthetics industry - where science, sustainability, and service go hand-in-hand.

with love
Di & Family

Technology meets Anatomy

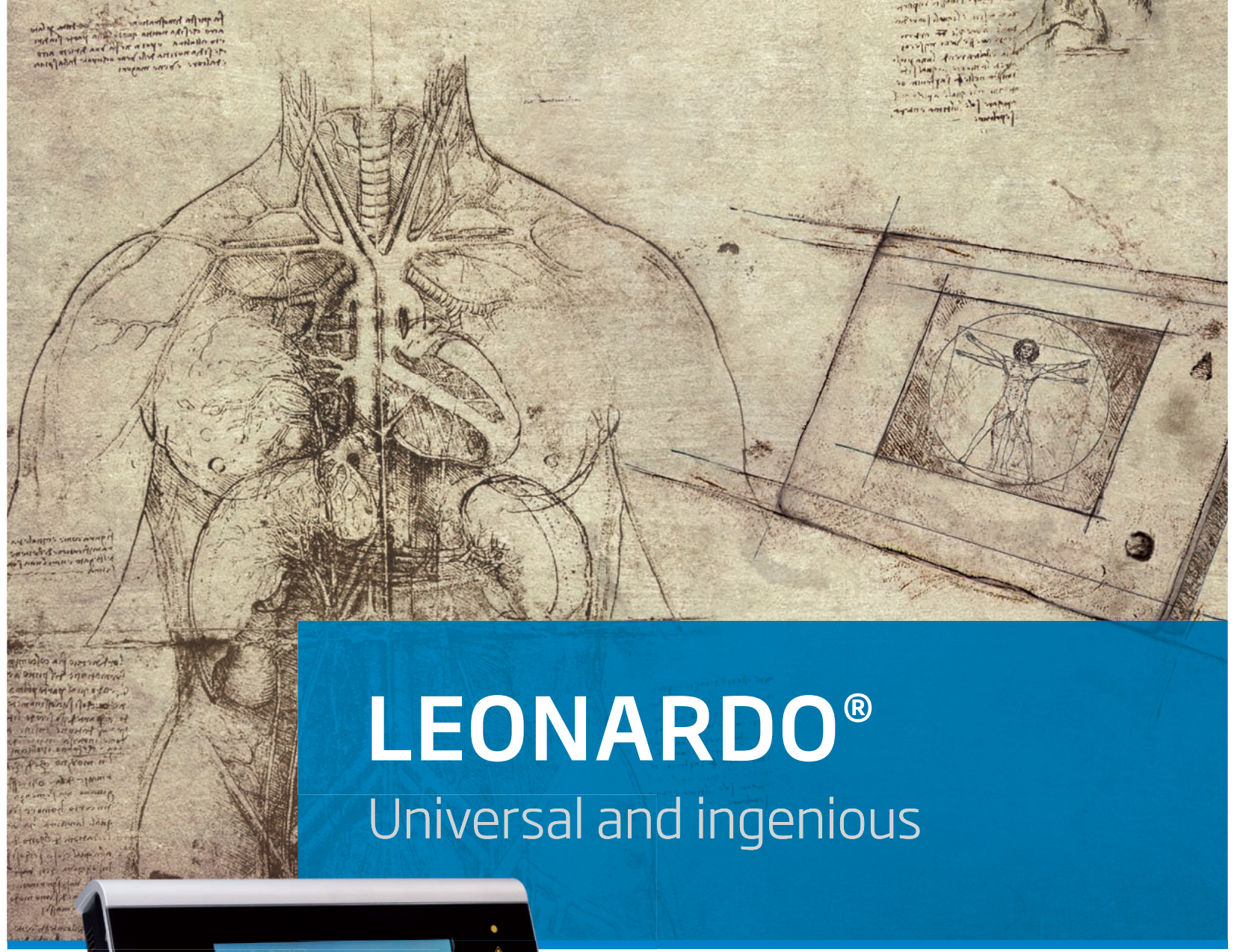

LEONARDO[®]
Universal and ingenious




- Cutting
- Vaporization
- Tissue shrinking
- Coagulation
- Dual wavelength laser
individually selected
or blended




A whole new world of therapeutic applications and clinical results

biolitec® presents the LEONARDO® Laser - the most versatile and universal medical laser in the market today. This highly compact diode laser features the combination of two wavelengths, 980 nm and 1470 nm offering a variety of tissue interactions. Each wavelength can be individually selected or **blended** together to offer the perfect desired tissue effects such as incision, excision, vaporization, hemostasis and coagulation of soft tissue with contact or non-contact delivery options for open and endoscopic procedures. For the first time the clinicians can perform laser surgery selectively, with settings individually tailored to the tissue type and the desired tissue effects and thus corresponding to the therapeutic needs.

The ability to choose a wavelength mix opens a whole new world of therapeutic applications and improves both the treatment outcomes for the patients and **extends** the clinicians experience and expertise. The LEONARDO® Laser is designed to work in perfect combination with a broad spectrum of special medical fibers and application kits developed by biolitec® and its companies. biolitec® is one of the **world's most** technologically advanced suppliers of fiberoptic products. The biolitec® treatment methods are routinely performed and validated in highly respected medical centers worldwide and are the number one choice for treating a wide variety of diseases and medical conditions.

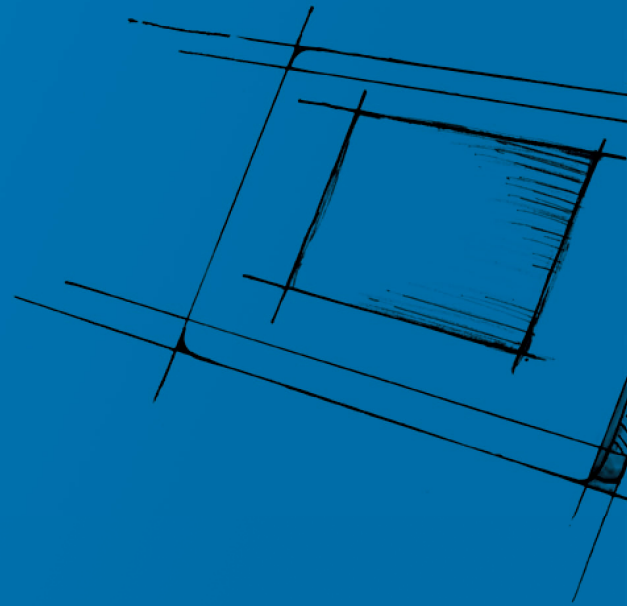

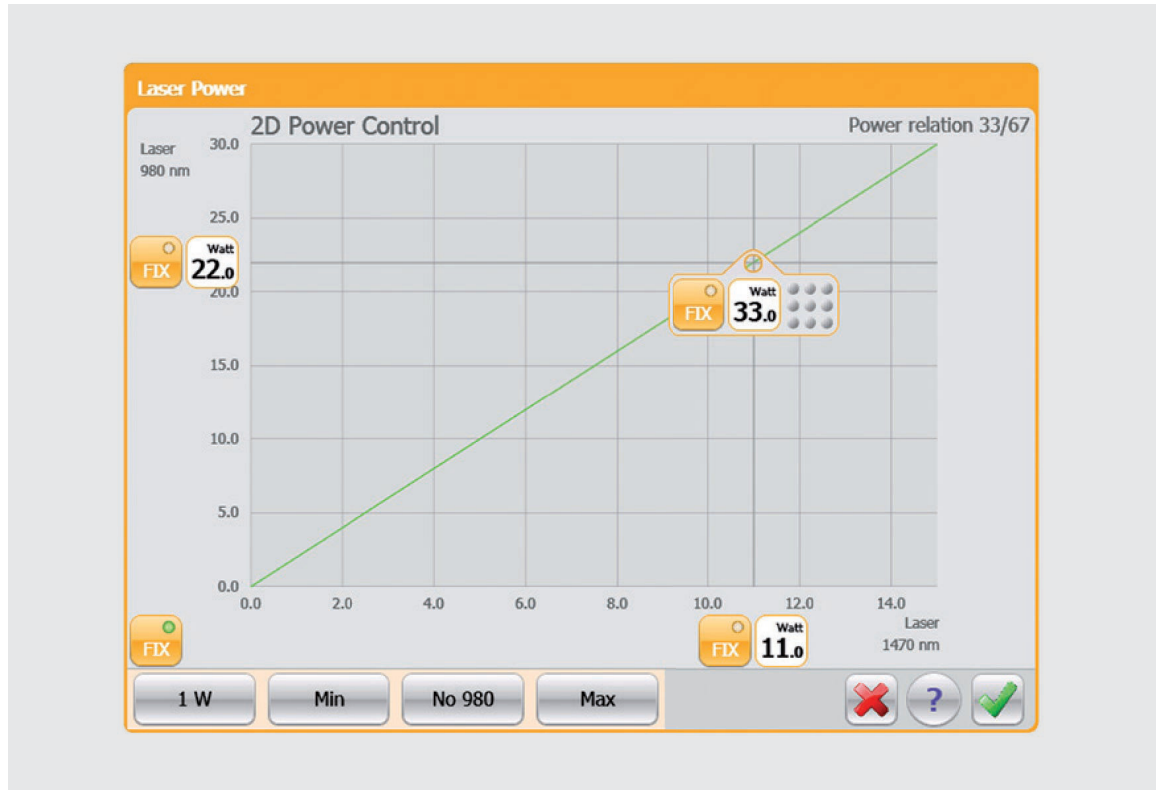

2D Power Control

LEONARDO®'s intuitive 2D Power Control enables the user to choose a combination of different wavelengths and power settings with a simple touch of the screen.

New Fiber Connector

The new fiber connector facilitates to plug the fiber into the laser. It is equipped with an electronic signature for increased patient safety. It prevents usage beyond the product's lifetime and hazards caused by inserting unsuitable fibers to the laser.

Advantages

Versatile and universal


- Broad spectrum of minimally invasive therapeutic laser applications

User-Friendly

- Intuitive use with touch screen and fast set-up
- Selection between pre-set modes or individualized settings
- Choice between green or red aiming beam

Economic

- Two wavelength in one compact and space-saving laser system
- Multidisciplinary use
- Low-maintenance and reliable laser diodes

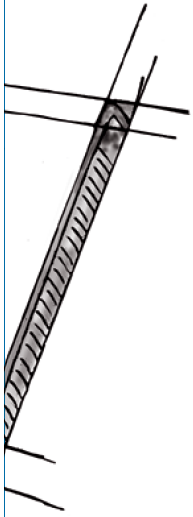

LipoLas[®] – Laser Lift

Minimally invasive laser therapy for skin tightening and body contouring



ONE PLANET
AESTHETICS

- Significant fat reduction and tissue tightening
- Precise and gentle fiber guidance
- Outpatient treatment
- Versatile in use

Applications for fat reduction and tissue tightening




FACE




CHIN




NECK

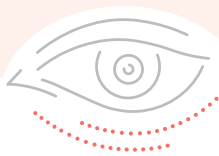


EYE

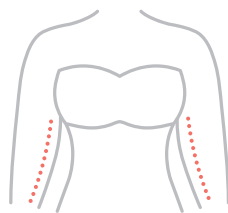


UPPER ARMS

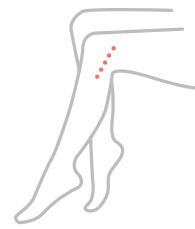


KNEE

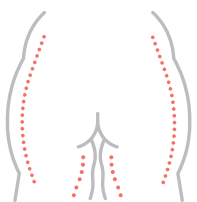

BOTTOM

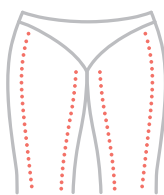

THIGHS

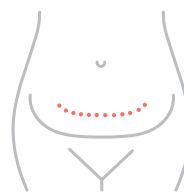

BELLY

What is LipoLas[®] – Laser Lift?


LipoLas[®] - Laser Lift is a minimally invasive procedure for skin tightening and fat reduction without scarring. The procedure can be performed on an outpatient basis under local anesthesia, e.g. TLA, or in combination with liposuction. A thin laser fiber is inserted subcutaneously into the tissue to be treated. No skin incision is required for this. The desired area is treated evenly by slow and fan-shaped release of laser energy.

The position of the laser fiber is continuously controlled with the pilot beam visible through the skin. The specific interaction of the laser light in the tissue ensures a controlled, precise application and impresses with good results and long-term effects.

- Firms and tightens the skin
- Reduces excess fat
- Defines contours on the face and other areas of the body

Mode of action

The LEONARDO[®] is optionally equipped with 2 wavelengths 1470nm and 980nm, which have a high absorption in water, fat and hemoglobin. The complementary and simultaneous use of both wavelengths enables optimum reduction of excess fatty tissue and excellent tightening with maximum protection of the surrounding tissue. Due to the laser-induced sclerosis of the small blood vessels and stimulation of the formation of new collagen fibers, the postoperative course is also favored.

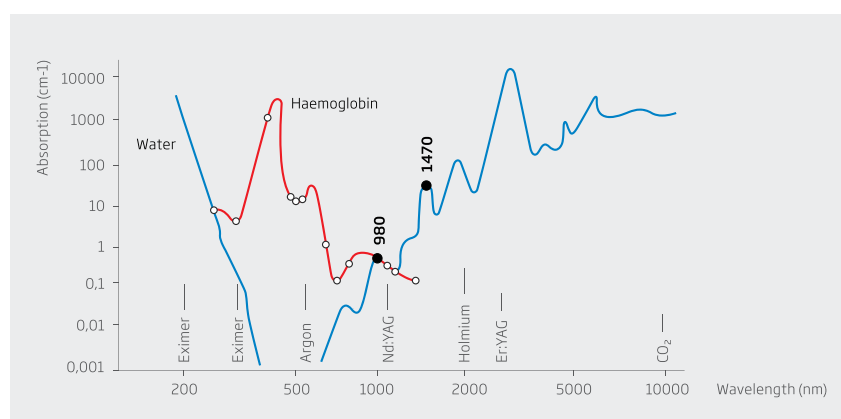

Laser absorption in hemoglobin and water

atraumatic fiber tip



Expansion of the therapeutic spectrum

Whether in the operating room of a large hospital, as an outpatient in a surgical center or in a private practice - the use of the LEONARDO® laser system significantly expands the therapeutic spectrum for the physician.

Due to the described mode of action, almost all areas of the body can be tightened and defined with LipoLas®, regardless of whether they are small and sensitive areas around the eyes, face and neck or large areas on the body. LipoLas® is also ideal for the treatment of dents in the skin, known as cellulite.

In addition to purely aesthetic applications, the LEONARDO® laser system also offers other possible applications, such as in phlebology, proctology, gynecology (LaEvita*) and the treatment of hidradenitis suppurativa.

* LaEvita is the new, gentle solution for the treatment of vaginal atrophy and mild stress incontinence.

atraumatic
fiber tip

LipoLas® -
Laser Lift
Fiber

Advantages of LipoLas® - Laser Lift

- Minimally invasive procedure without scarring
- Effective skin tightening with long-term effect
- Removal of excess fat from small areas of skin
- Definition of face and body areas
- Can be used on almost all skin areas



LaEvita

Regain your femininity

The laser treatment of genito-urinary syndrome of the menopause

biolitec® has the complete solution to the main disorders affecting the female reproductive system.

Why LaEvita?

LaEvita is the new virtually pain-free solution to the treatment of **vaginal atrophy** and **mild stress urinary incontinence**. Unlike other CO₂-lasers, this is a **NON-ablative and NON-destructive tissue technique** that is virtually painless and well tolerated by patients.

How does it work?

Thanks to the special applicator, the LEONARDO® Dual laser emits a gentle, painless light that stimulates collagen production and restores vaginal elasticity and hydration, and improves vascularization of the vaginal mucosa. Each session lasts about 10 minutes and is performed without anesthesia. The treatment includes 3 sessions at 4-week intervals. In more severe cases, 1 additional session may be necessary. The procedure is performed on an outpatient basis, is virtually painless for the woman, and does not require hospitalization.

What is required?

The technique requires the use of the LEONARDO® Dual laser, the use of special reusable LaEvita handpieces and dedicated kits for the treatment of vaginal atrophy or urinary incontinence

Advantages

- It significantly reduces vaginal dryness and pain without side effects
- It helps maintain a fulfilling intimate life
- Fully outpatient procedure
- No anesthesia required
- No side effects
- Effective
- Almost no pain
- No downtime



LaEvita Handpieces*

REF	Product	Outer ø [mm]	Inner ø [mm]	Length ø [mm]
400500400	LaEvita Small	24	2,20	100
400500410	LaEvita Medium	24	2,20	135
400500420	LaEvita Large	28	2,20	135
400500430	LaEvita Handpieces Set - (LaEvita small, medium, large)			

*The Vaginal Handpiece is a reusable for use in gynecological vaginal treatments.



LEONARDO® DUAL 45
The complete solution
for gynecology and more



Kits

REF	Product	Components	Quantity
503200610	Vaginal Atrophy LaEvita Kit, IC	Coupling Fiber, IC	1
		Side Fiber Tip	10
503200620	Urinary Incontinence LaEvita Kit, IC	Coupling Fiber, IC	1
		Radial Fiber Tip	10
		Side Fiber Tip	10



ONE PLANET
AESTHETICS

FROM
Our Family
TO YOURS

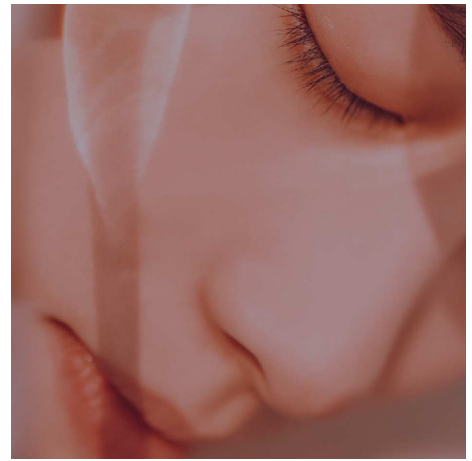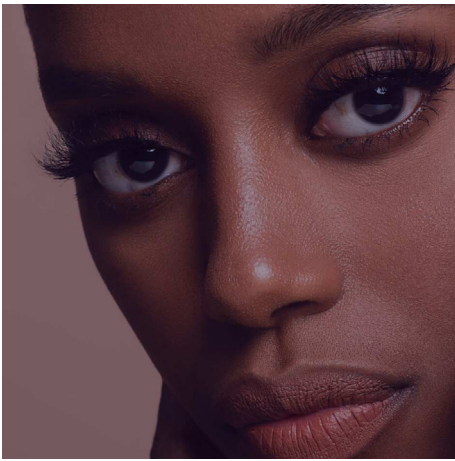

☎ 0800 772 0402

🌐 www.oneplanetaesthetics.co.uk

✉ info@oneplanetaesthetics.co.uk

📷 📘 @oneplanetaesthetics



St Clare House, Ninth Floor, Greyfriars,
Princes Street, Ipswich, Suffolk, IP1 1LP